H2401

Maia Gal

Galvanic / Corriente galvánica / Courant galvanique
Corrente galvanica / Galvanischer Strom

EN	2
ES	5
FR	9
IT	13
DE	17

DDUUEETT EUROPE would like to thank you for choosing Sonia products. We hope that this instruction manual will help you get the most out of this product.

If you have any doubts regarding its use and/or operation, please contact your distributor.

BEFORE USE

- Please read the entire instruction manual carefully before using the product. Take note of any warnings included in this manual.
- Always keep this instruction manual handy and if you give the device to another person, hand over the manual as it is an essential part of the product.
- Some components of the product may differ from the image that appears in catalogues, websites, brochures, etc.
- In order to enhance the product, the manufacturer reserves the right to modify its features and/or components without prior notice.

CORRECT USE OF THE DEVICE

- This product has been designed to be used by beauty professionals at beauty salons and should not be used for other purposes.
- Please note that any misuse of the device entails the loss of the guarantee.
- Do not make any modification to the device without the authorisation of DDUUEETT EUROPE SL.
- DDUUEETT EUROPE SL shall not be held liable for any material or personal damage due to the misuse of the product.

CONTRAINDICATIONS

- This product should not be used on:
 - People with metal implants in their body.
 - People with pacemakers or heart problems.
 - Pregnant women.
 - People with sensitive skin, lesions, inflammations, or any skin problem or disease.
- With people undergoing medical treatment or a diagnostic process, the product must be used with medical authorisation.

PRECAUTIONS

- Before using it for the first time and after each use, check the device and the mains cable for any damage.
- Unplug the device from the mains before starting to clean it or undertaking any type of maintenance.
- Do not use the product outdoors or in places

where it may be affected by humidity or inclement weather.

- If the power supply cable is defective, unplug the appliance immediately and inform your dealer or specialised technician.
- Position the power supply cable properly so that people cannot trip over it and make sure you plug it into an earthed socket.
- At the end of the working day, disconnect the device and unplug the power supply cable from the mains.
- If you detect any anomaly in the operation of the product, please contact your distributor or specialised technician.

SAFETY INDICATIONS

- The device and accessories should be kept out of the reach of children.
- Only plug the device into a properly installed socket. The local mains voltage must correspond to the technical specifications of the device.
- The power socket must be accessible in case you have to quickly unplug the device at any given time.
- Always pull out the actual plug and not the mains cable.
- Never tamper with any electrical component or open up/repair any component of the device.
- Only your distributor or specialised technician should repair the device.
- Do not switch on the device if there is visible damage to it or to the mains cable.
- Never touch the device or the mains cable with wet hands.


CONTENTS

- Main body x 1 unit
- Power cable x 1 unit
- Handpiece with double stick probe x 1 unit
- Handpiece with roller probe x 1 unit
- Instruction manual x 1 unit

TECHNICAL SPECIFICATIONS


Product	Galvanic
Reference	H2401
Voltage	85-265 V
Input frequency	50/60 Hz
Power	20 W
Dimensions	28 x 17 x 11 cm
Net weight	1,95 Kg

COMPONENT IDENTIFICATION

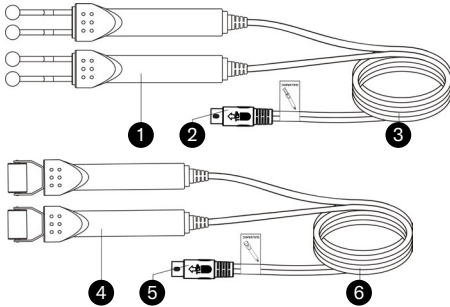- 1 Handpiece with double stick probe
 - 2 Double stick probe cable connector
 - 3 Double stick probe flexible cable
 - 4 Handpiece with roller probe
 - 5 Roller probe handpiece cable connector
 - 6 Roller probe flexible cable
- Double stick probe: to treat small or irregular areas of the skin
 - Roller probe: to treat large and regular areas of the skin


CONTROL PANEL

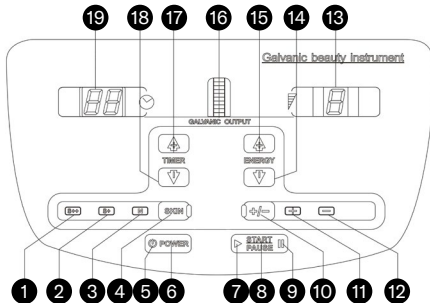- 1 Very sensitive skin indicator.
- 2 Sensitive skin indicator
- 3 Normal skin indicator
- 4 Skin type selection
- 5 Power on indicator
- 6 On / off button
- 7 Run indicator
- 8 Start / Stop button
- 9 Pause indicator

- 10 Ion conduction selection button
- 11 Cation conduction indicator (positive charge)
- 12 Anion conduction indicator (negative charge)
- 13 Power display
- 14 Decrease power button
- 15 Increase power button
- 16 Galvanic current conduction indicator
- 17 Increase time button
- 18 Decrease time button
- 19 Timer


REAR PANEL

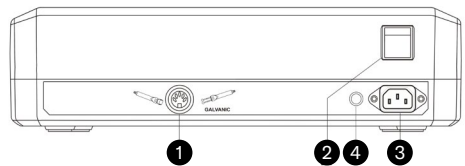- 1 Probe handpiece cable connection
- 2 Power supply switch
- 3 Electrical power supply connector
- 4 Fuse

CONNECTION INSTRUCTIONS

- Connect the power cable to the rear panel of the main body.
- Connect the cable of the handpiece to the corresponding output on the rear panel.
- It is only possible to connect one of the two handpieces at the same time, so you must select the one you need according to the treatment to be performed.

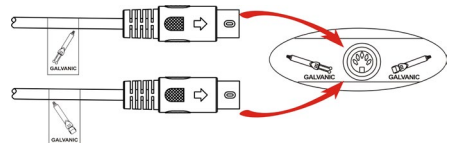POWER AND CONFIGURATION

- Turn on the power switch on the rear of the machine and press **POWER** on the control panel. The power indicator will light up. The machine will be in standby mode and ready for operation. On the control panel:
 - The timer will show 15 minutes.

- The power indicator will show 1.
- The normal skin indicator  will light up.
- The cation conduction (positive charge) indicator  will light up.
- On the  button, the pause indicator  will light up.
- Set the timer using the  or  buttons to choose the operating time of the machine, in a range of between 1 and 30 minutes.
- Set the power using the  or  buttons to choose the power of the machine, in a range of between 1 and 9.
- Skin type selection: use the  button to select the type of skin to be treated.
 - When the indicator  lights up, the selected skin type is normal.
 - When the  indicator lights up, the selected skin type is sensitive.
 - When the  indicator lights up, the selected skin type is very sensitive.
- Ion conduction selection: press  to select the polarisation of the probe with a red cable. The probe with a black cable will always have the polarity inverse to the one selected for the red cable.
 - One feature of this device is that the patient should not hold any electrodes to close the circuit, as one of the probes acts as positive and the other probe acts as negative, thus closing the circuit.
 - When the  indicator lights up, it indicates cation conduction (positive charge).
 - When the  indicator lights up, it indicates anion conduction (negative charge).
- Use the  button to start or stop the device.
 - When running, the  indicator will light up.
 - When paused or on standby, the  indicator will light up.

OPERATION

This device can be used with the roller probe handpieces for large or regular areas of the skin, or with the double stick probes for small or irregular areas of the skin.

- Connect the device and leave it in standby mode with the  indicator lit.
- Apply conductive gel or cream after cleaning the skin to be treated.
- Select the type of skin to be treated, the operating time and set the power to minimum.
- In order to correctly select the polarisation of the ion conduction with the  button, the pH value of the previously applied product should be considered, selecting one or the other polarity depending on whether it is an alkaline or an acidic

product.

- Press  and the run indicator  will light up.
- Take the two probes (or the roller or double stick ones) without letting them touch.
- Start to move them slowly and gradually set the power until you reach the level necessary for the treatment, always maintaining a slight stinging sensation.

When applying this treatment on the face, it is advisable to alternate the two probes (positive and negative) on both sides of the face in order to make equal effects on both sides of the face

Bear in mind:

- When the polarity is set for positive electrodes, there will be a stinging sensation given off by the roller probe that is stronger than the one given off by the negative electrode. Select it according to your needs.
- The selection of ion conduction to be used should coincide with the pH of the product applied.
 - When the probe with a red cable acts as a positive charge, the conductive solution should be acidic (+).
 - When the probe with a red cable acts as a negative charge, the conductive solution should be alkaline (-).
 - If the pH is not indicated in the solution to be used, use a pH indicator strip to determine it.
- If the run indicator  is flashing, the circuit is not closed and there is no ion conduction.

CLEANING AND MAINTENANCE

- Before cleaning the device, unplug it from the mains.
- Clean the device with a soft, dry cloth. Do not use liquids to clean it.
- After each use, the various components and accessories must be sterilised and left to dry.
- For maintenance or repairs of any part of the product, please contact your distributor or specialised technician.

RECYCLING

If you need to dispose of the product or any of its parts or accessories, please check and respect the laws and regulations applicable to waste management and environmental protection.

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